

## GMP FITC anti-human CD66b Antibody

<b>Catalog# / Size</b>	260188 / 100 tests
<b>Clone</b>	G10F5
<b>Workshop</b>	VI MA81
<b>Other Names</b>	CD67, CGM6, NCA-95, CEACAM8
<b>Isotype</b>	Mouse IgM, κ
<b>Description</b>	CD66b is a 95-100 kD glycosylphosphatidylinositol (GPI)-linked protein also known as CD67, CGM6, and NCA-95. CD66b is a member of the immunoglobulin superfamily, carcinoembryonic antigen (CEA)-like subfamily. CD66b, expressed on granulocytes, has been reported to induce activation in neutrophils and to be involved in heterophilic adhesion with CD66c.

### Product Details

<b>Reactivity</b>	Human
<b>Antibody Type</b>	Monoclonal
<b>Host Species</b>	Mouse
<b>Formulation</b>	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and 0.2% (w/v) BSA (origin USA) and a stabilizer.
<b>Preparation</b>	The antibody was purified by affinity chromatography and conjugated with FITC under optimal conditions.
<b>Concentration</b>	200 µg/mL
<b>Storage &amp; Handling</b>	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. <b>Do not freeze.</b>
<b>Application</b>	<a href="#">FC - Quality tested</a>
<b>Recommended Usage</b>	Each lot of this antibody is quality control tested by <a href="#">immunofluorescent staining with flow cytometric analysis</a> . For flow cytometric staining, the suggested use of this reagent is 5 µL per million cells in 100 µL staining volume or 5 µL per 100 µL of whole blood. It is recommended that the reagent be titrated for optimal performance for each application.
<b>Excitation Laser</b>	Blue Laser (488 nm)
<b>Application Notes</b>	Additional reported applications (for the relevant formats) include: immunohistochemical staining of acetone-fixed frozen, formalin-fixed paraffin-embedded tissue sections, and spatial biology (IBEX) <sup>5,6</sup> .
<b>Application References</b>	<ol style="list-style-type: none"> <li>Schlossman S, <i>et al.</i> Eds. 1995. Leucocyte Typing V. Oxford University Press. New York.</li> <li>Kishimoto T, <i>et al.</i> Eds. 1997. Leucocyte Typing VI. Garland Publishing Inc. London.</li> <li>Norling LV, <i>et al.</i> 2012. <i>Arterioscler Thromb Vasc Biol.</i> 32:1970. <a href="#">PubMed</a></li> <li>Meinke P, <i>et al.</i> 2015. <i>Neuroimmunol Discord.</i> 25:127. <a href="#">PubMed</a></li> <li>Radtke AJ, <i>et al.</i> 2020. <i>Proc Natl Acad Sci USA.</i> 117:33455-33465. (SB) <a href="#">PubMed</a></li> <li>Radtke AJ, <i>et al.</i> 2022. <i>Nat Protoc.</i> 17:378-401. (SB) <a href="#">PubMed</a></li> </ol>
<b>(PubMed link indicates BioLegend citation)</b>	
<b>Disclaimer</b>	<p><b>GMP RUO Flow Cytometry Antibodies.</b> BioLegend GMP RUO fluorophore conjugated antibodies are manufactured in a dedicated GMP facility and compliant with ISO 13485:2016. For research use only. Not for use in diagnostic or therapeutic procedures. Our processes include:</p> <ul style="list-style-type: none"> <li>• Batch-to-batch consistency</li> <li>• Material traceability</li> <li>• Documented procedures</li> <li>• Documented employee training</li> <li>• Equipment maintenance and monitoring records</li> <li>• Lot-specific certificates of analysis</li> <li>• Quality audits per ISO 13485:2016</li> </ul>

- QA review of released products

## Antigen Details

---

<b>Structure</b>	Ig superfamily, CEA antigen group, GPI-linked glycoprotein, 95-100 kD
<b>Distribution</b>	Granulocytes
<b>Function</b>	Cell adhesion, neutrophil activation
<b>Ligand/Receptor</b>	CD66c
<b>Cell Type</b>	Granulocytes, Neutrophils
<b>Biology Area</b>	Immunology
<b>Molecular Family</b>	Adhesion Molecules, CD Molecules
<b>Antigen References</b>	1. Kuijpers T, <i>et al.</i> 1993. <i>J. Immunol.</i> 151:4934. 2. Kuroki M, <i>et al.</i> 1992. <i>J. Leuk. Biol.</i> 52:551.
<b>Gene ID</b>	<a href="#">1088</a>

## Related Protocols

---

[Cell Surface Flow Cytometry Staining Protocol](#)


## Other Formats

---

FITC anti-human CD66b, Purified anti-human CD66b, Pacific Blue™ anti-human CD66b, PE anti-human CD66b, PerCP/Cyanine5.5 anti-human CD66b, Alexa Fluor® 647 anti-human CD66b, Alexa Fluor® 700 anti-human CD66b, PE/Cyanine7 anti-human CD66b, APC anti-human CD66b, Biotin anti-human CD66b, PE/Dazzle™ 594 anti-human CD66b, Alexa Fluor® 594 anti-human CD66b, APC/Cyanine7 anti-human CD66b

## Product Data

---

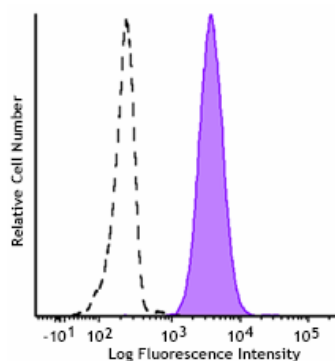

Typical results from human peripheral blood granulocytes stained either with G10F5 FITC used at 5  $\mu$ L/test (filled histogram) or with an isotype control (open histogram).

For research use only. Not for resale. BioLegend will not be held responsible for patent infringement or other violations that may occur with the use of our products.

\*These products may be covered by one or more Limited Use Label Licenses (see the BioLegend Catalog or our website, [www.biolegend.com/ordering#license](http://www.biolegend.com/ordering#license)). BioLegend products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products, reverse engineer functionally similar materials, or to provide a service to third parties without written approval of BioLegend. By use of these products you accept the terms and conditions of all applicable Limited Use Label Licenses.

BioLegend Inc., 8999 BioLegend Way, San Diego, CA 92121 [www.biolegend.com](http://www.biolegend.com)  
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587