

Brilliant Violet 605™ anti-human CD63 Antibody

Catalog# / Size	353047 / 25 tests 353048 / 100 tests
Clone	H5C6
Regulatory Status	RUO
Workshop	HCDM listed
Other Names	LIMP, LAMP-3, Melanoma-associated antigen (ME491), PItgp40, gp55
Isotype	Mouse IgG1, κ
Description	CD63 is a 53 kD type III lysosomal glycoprotein also known as LIMP, LAMP-3, gp55, and melanoma-associated antigen (ME491). CD63 is a member of the tetraspan transmembrane superfamily (TM4SF) protein expressed on activated platelets, monocytes/macrophages, endothelium, fibroblasts, osteoclasts, and smooth muscle cells. CD63 may be involved in platelet activation and is thought to function as a transmembrane adaptor protein. CD63 has been shown to associate with CD9, CD81, VLA-3, and VLA-6.

Product Details

Verified Reactivity	Human, African Green, Baboon, Cynomolgus, Rhesus
Antibody Type	Monoclonal
Host Species	Mouse
Immunogen	T cell line HPB-ALL
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA)
Preparation	The antibody was purified by affinity chromatography and conjugated with Brilliant Violet 605™ under optimal conditions.
Concentration	Lot-specific (to obtain lot-specific concentration, please enter the lot number in our Concentration and Expiration Lookup or Certificate of Analysis online tools.)
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	<p>Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis. For flow cytometric staining, the suggested use of this reagent is 5 µL per million cells in 100 µL staining volume or 5 µL per 100 µL of whole blood. It is recommended that the reagent be titrated for optimal performance for each application.</p> <p>Brilliant Violet 605™ excites at 405 nm and emits at 603 nm. The bandpass filter 610/20 nm is recommended for detection, although filter optimization may be required depending on other fluorophores used. Be sure to verify that your cytometer configuration and software setup are appropriate for detecting this channel. Refer to your instrument manual or manufacturer for support. Brilliant Violet 605™ is a trademark of Sirigen Group Ltd.</p> <p>Learn more about Brilliant Violet™.</p> <p>This product is subject to proprietary rights of Sirigen Inc. and is made and sold under license from Sirigen Inc. The purchase of this product conveys to the buyer a non-transferable right to use the purchased product for research purposes only. This product may not be resold or incorporated in any manner into another product for resale. Any use for therapeutics or diagnostics is strictly prohibited. This product is covered by U.S. Patent(s), pending patent applications and foreign equivalents.</p>
Excitation Laser	Violet Laser (405 nm)
Application Notes	Additional reported applications (for the relevant formats) include: Western blotting ¹ ,

immunofluorescence², and immunoprecipitation¹.

Application References

(PubMed link indicates BioLegend citation)

1. Hildreth JE, *et al.* 1991. *Blood* 77:121. (IP, WB)
2. Beatty WL, *et al.* 2006. *J. Cell Sci.* 119:350. (IF)

RRID

AB_2890780 (BioLegend Cat. No. 353047)
AB_2888782 (BioLegend Cat. No. 353048)

Antigen Details

Structure	Tetraspan transmembrane superfamily (TM4SF), type III lysosomal glycoprotein, 53 kD
Distribution	Activated platelets, monocytes, macrophages, endothelium, fibroblasts, osteoclasts, and smooth muscle
Function	Platelet activation
Interaction	CD9, CD81, VLA-3, and VLA-6
Cell Type	B cells, Endothelial cells, Fibroblasts, Macrophages, Monocytes, Osteoclasts, Platelets
Biology Area	Cell Adhesion, Cell Biology, Immunology
Molecular Family	Adhesion Molecules, CD Molecules

Antigen References

1. Azorsa DO, *et al.* 1991. *Blood* 78:280.
2. Kishimoto T, *et al.* Eds. 1997. *Leukocyte Typing V1*. Oxford University Press New York.
3. Hildreth JE, *et al.* 1991. *Blood* 77:121.
4. Anzai N, *et al.* 2002. *Blood* 99:4413.

Gene ID

[967](#)


Related Protocols

[Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

PE anti-human CD63, FITC anti-human CD63, Pacific Blue™ anti-human CD63, APC anti-human CD63, PE/Cyanine7 anti-human CD63, Alexa Fluor® 647 anti-human CD63, Biotin anti-human CD63, PerCP/Cyanine5.5 anti-human CD63, PE/Dazzle™ 594 anti-human CD63, Brilliant Violet 421™ anti-human CD63, Brilliant Violet 650™ anti-human CD63, APC/Fire™ 750 anti-human CD63, Alexa Fluor® 700 anti-human CD63, Brilliant Violet 510™ anti-human CD63, Alexa Fluor® 594 anti-human CD63, TotalSeq™-A0404 anti-human CD63, Alexa Fluor® 488 anti-human CD63, Ultra-LEAF™ Purified anti-human CD63, Brilliant Violet 711™ anti-human CD63, Brilliant Violet 785™ anti-human CD63, APC/Cyanine7 anti-human CD63, TotalSeq™-B0404 anti-human CD63, TotalSeq™-C0404 anti-human CD63, Brilliant Violet 605™ anti-human CD63 Antibody, PE/Cyanine5 anti-human CD63 Antibody

Product Data

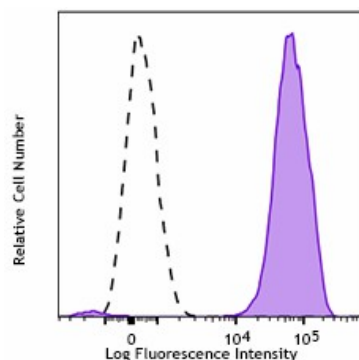

Thrombin-activated human peripheral blood platelets were stained with CD63 (clone H5C6) Brilliant Violet 605™ (filled histogram) or mouse IgG1, κ Brilliant Violet 605™ isotype control (open histogram).

For research use only. Not for diagnostic use. Not for resale. BioLegend will not be held responsible for patent infringement or other violations that may occur with the use of our products.

*These products may be covered by one or more Limited Use Label Licenses (see the BioLegend Catalog or our website, www.biolegend.com/ordering#license). BioLegend products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products, reverse engineer functionally similar materials, or to provide a service to third parties without written approval of BioLegend. By use of these products you accept the terms and conditions of all applicable Limited Use Label Licenses. Unless otherwise indicated, these products are for research use only and are not intended for human or animal diagnostic, therapeutic or commercial use.

BioLegend Inc., 8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587