

## Brilliant Violet 785™ anti-human CD366 (Tim-3) Antibody

<b>Catalog# / Size</b>	345031 / 25 tests 345032 / 100 tests
<b>Clone</b>	F38-2E2
<b>Regulatory Status</b>	RUO
<b>Workshop</b>	HCDM listed
<b>Other Names</b>	T cell immunoglobulin and mucin domain containing protein 3, hepatitis virus cellular receptor 2, CD366
<b>Isotype</b>	Mouse IgG1, κ
<b>Description</b>	CD366 (Tim-3) is a transmembrane protein also known as T cell immunoglobulin and mucin domain containing protein-3. Tim-3 is expressed at high levels on activated T cells (preferentially on Th1 cells, monocytes/macrophages, and dendritic cells). Tim-3 has also been shown to exist as a soluble protein. Cells expressing Tim-3 are present at high levels in the CNS of animals at the onset of experimental autoimmune encephalomyelitis (EAE), a disease mediated by lymphocytes secreting Th1-like cytokines. Tim-3 has been proposed to inhibit Th1-mediated immune responses and promote immunological tolerance.

### Product Details

---

<b>Verified Reactivity</b>	Human
<b>Antibody Type</b>	Monoclonal
<b>Host Species</b>	Mouse
<b>Immunogen</b>	Human Tim-3 fusion protein
<b>Formulation</b>	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA).
<b>Preparation</b>	The antibody was purified by affinity chromatography and conjugated with Brilliant Violet 785™ under optimal conditions.
<b>Concentration</b>	Lot-specific (to obtain lot-specific concentration, please enter the lot number in our <a href="#">Concentration and Expiration Lookup</a> or <a href="#">Certificate of Analysis</a> online tools.)
<b>Storage &amp; Handling</b>	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. <b>Do not freeze.</b>
<b>Application</b>	<a href="#">FC - Quality tested</a>
<b>Recommended Usage</b>	Each lot of this antibody is quality control tested by <a href="#">immunofluorescent staining with flow cytometric analysis</a> . For flow cytometric staining, the suggested use of this reagent is 5 µl per million cells in 100 µl staining volume or 5 µl per 100 µl of whole blood.

Brilliant Violet 785™ excites at 405 nm and emits at 785 nm. The bandpass filter 780/60 nm is recommended for detection, although filter optimization may be required depending on other fluorophores used. **Be sure to verify that your cytometer configuration and software setup are appropriate for detecting this channel.** Refer to your instrument manual or manufacturer for support. Brilliant Violet 785™ is a trademark of Sirigen Group Ltd.

[Learn more about Brilliant Violet™.](#)

This product is subject to proprietary rights of Sirigen Inc. and is made and sold under license from Sirigen Inc. The purchase of this product conveys to the buyer a non-transferable right to use the purchased product for research purposes only. This product may not be resold or incorporated in any manner into another product for resale. Any use for therapeutics or diagnostics is strictly prohibited. This product is covered by U.S. Patent(s), pending patent applications and foreign equivalents.

<b>Excitation Laser</b>	Violet Laser (405 nm)
-------------------------	-----------------------

<b>Application Notes</b>	Additional reported applications (for relevant formats of this clone) include: costimulation <sup>1</sup> (clone 2E2 has been shown to enhance T-cell receptor mediated activation and cytokine secretion) and blocking <sup>2,3</sup> .
<b>Application References</b>	<ol style="list-style-type: none"> <li>Hastings WD, <i>et al.</i> 2009. <i>Eur. J. Immunol.</i> 39:2492. (Costim)</li> <li>Jones RB, <i>et al.</i> 2008. <i>J. Exp. Med.</i> 205:2763. (Block)</li> <li>Klibi J, <i>et al.</i> 2009. <i>Blood</i> 113:1957. (FC, Block)</li> </ol>
<b>(PubMed link indicates BioLegend citation)</b>	
<b>Product Citations</b>	<ol style="list-style-type: none"> <li>Eyquem J, <i>et al.</i> 2017. <i>Nature.</i> 543:113-117. <a href="#">PubMed</a></li> <li>Amor C, <i>et al.</i> 2020. <i>Nature.</i> 583:127. <a href="#">PubMed</a></li> <li>Gibellini L, <i>et al.</i> 2020. <i>EMBO Mol Med.</i> 12:e13001. <a href="#">PubMed</a></li> <li>Kim Y, <i>et al.</i> 2021. <i>Oncol Lett.</i> 1.022222222. <a href="#">PubMed</a></li> <li>Mathewson ND, <i>et al.</i> 2021. <i>Cell.</i> 184(5):1281-1298.e26. <a href="#">PubMed</a></li> <li>Chow MT <i>et al.</i> 2019. <i>Immunity.</i> 50(6):1498-1512. <a href="#">PubMed</a></li> <li>Kachеровsky N, <i>et al.</i> 2019. <i>Nat Biomed Eng.</i> 0.66875. <a href="#">PubMed</a></li> <li>Karlsson J, <i>et al.</i> 2020. <i>Nat Commun.</i> 1.773611111. <a href="#">PubMed</a></li> </ol>
<b>RRID</b>	AB_2565832 (BioLegend Cat. No. 345031) AB_2565833 (BioLegend Cat. No. 345032)

## Antigen Details

<b>Structure</b>	Transmembrane protein containing immunoglobulin domain and mucin-like domain; can exist as a soluble form lacking mucin and transmembrane domains
<b>Distribution</b>	Activated T cells, preferentially on Th1 cells, monocytes, dendritic cells
<b>Function</b>	Plays a role in regulating macrophage activation, T cell apoptosis and immune tolerance
<b>Ligand/Receptor</b>	Galectin-9
<b>Cell Type</b>	Dendritic cells, Monocytes, T cells, Th1, Tregs
<b>Biology Area</b>	Immunology, Inhibitory Molecules
<b>Molecular Family</b>	CD Molecules, Immune Checkpoint Receptors
<b>Antigen References</b>	<ol style="list-style-type: none"> <li>Hafler DA and Kuchroo V. 2008. <i>J. Exp. Med.</i> 205:2699.</li> <li>Zhu C, <i>et al.</i> 2005. <i>Nat. Immunol.</i> 6:1245.</li> <li>Wang F, <i>et al.</i> 2009. <i>Immunobiology</i> 214:342.</li> </ol>
<b>Gene ID</b>	<a href="#">84868</a>


## Related Protocols

[Cell Surface Flow Cytometry Staining Protocol](#)

## Other Formats

Purified anti-human CD366 (Tim-3), PE anti-human CD366 (Tim-3), Brilliant Violet 421™ anti-human CD366 (Tim-3), Ultra-LEAF™ Purified anti-human CD366 (Tim-3), APC anti-human CD366 (Tim-3), PE/Cyanine7 anti-human CD366 (Tim-3), PerCP/Cyanine5.5 anti-human CD366 (Tim-3), Brilliant Violet 605™ anti-human CD366 (Tim-3), FITC anti-human CD366 (Tim-3), Purified anti-human CD366 (Tim-3) (Maxpar® Ready), Brilliant Violet 711™ anti-human CD366 (Tim-3), APC/Cyanine7 anti-human CD366 (Tim-3), Brilliant Violet 785™ anti-human CD366 (Tim-3), Brilliant Violet 650™ anti-human CD366 (Tim-3), Brilliant Violet 510™ anti-human CD366 (Tim-3), PE/Dazzle™ 594 anti-human CD366 (Tim-3), GolnVivo™ Purified anti-human CD366 (Tim-3), APC/Fire™ 750 anti-human CD366 (Tim-3), Pacific Blue™ anti-human CD366 (Tim-3), Biotin anti-human CD366 (Tim-3), TotalSeq™-A0169 anti-human CD366 (Tim-3), TotalSeq™-C0169 anti-human CD366 (Tim-3), PE/Cyanine5 anti-human CD366 (Tim-3), TotalSeq™-B0169 anti-human CD366 (Tim-3), Brilliant Violet 750™ anti-human CD366 (Tim-3) Antibody, TotalSeq™-D0169 anti-human CD366 (Tim-3), PE/Fire™ 640 anti-human CD366 (Tim-3)

## Product Data

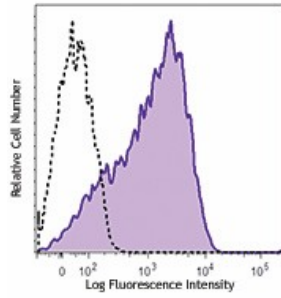

PHA-stimulated (three days) human peripheral blood lymphocytes were stained with anti-human CD366 (clone F38-2E2) Brilliant Violet 785™ (filled histogram) or mouse IgG1, κ Brilliant Violet 785™ (open histogram).

For research use only. Not for diagnostic use. Not for resale. BioLegend will not be held responsible for patent infringement or other violations that may occur with the use of our products.

\*These products may be covered by one or more Limited Use Label Licenses (see the BioLegend Catalog or our website, [www.biolegend.com/ordering#license](http://www.biolegend.com/ordering#license)). BioLegend products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products, reverse engineer functionally similar materials, or to provide a service to third parties without written approval of BioLegend. By use of these products you accept the terms and conditions of all applicable Limited Use Label Licenses. Unless otherwise indicated, these products are for research use only and are not intended for human or animal diagnostic, therapeutic or commercial use.

BioLegend Inc., 8999 BioLegend Way, San Diego, CA 92121 [www.biolegend.com](http://www.biolegend.com)  
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587