

Brilliant Violet 510™ anti-human CD134 (OX40) Antibody

Catalog# / Size	350025 / 25 tests 350026 / 100 tests
Clone	Ber-ACT35 (ACT35)
Regulatory Status	RUO
Other Names	ACT35 antigen, tumor necrosis factor receptor superfamily member 4 (TNFRSF4)
Isotype	Mouse IgG1, κ
Description	CD134, also known as OX40 and TNFRSF4, is a 50 kD type I transmembrane glycoprotein. It is a member of the TNF receptor family. OX40 is expressed on activated T lymphocytes including Th1, Th2, Th17, and Treg cells. The interaction of OX40 with OX40L results in B cell proliferation and antibody secretion, regulation of primary T cell expansion, and T cell survival. OX40 influences the size of the T cell memory pool and regulation of CD4 ⁺ T cell tolerance.

Product Details

Verified Reactivity	Human
Reported Reactivity	Chimpanzee
Antibody Type	Monoclonal
Host Species	Mouse
Immunogen	HTLV 1-transformed HUT 102 cells
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA).
Preparation	The antibody was purified by affinity chromatography and conjugated with Brilliant Violet 510™ under optimal conditions.
Concentration	Lot-specific (to obtain lot-specific concentration, please enter the lot number in our Concentration and Expiration Lookup or Certificate of Analysis online tools.)
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	<p>Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis. For flow cytometric staining, the suggested use of this reagent is 5 μl per million cells in 100 μl staining volume or 5 μl per 100 μl of whole blood.</p> <p>Brilliant Violet 510™ excites at 405 nm and emits at 510 nm. The bandpass filter 510/50 nm is recommended for detection, although filter optimization may be required depending on other fluorophores used. Be sure to verify that your cytometer configuration and software setup are appropriate for detecting this channel. Refer to your instrument manual or manufacturer for support. Brilliant Violet 510™ is a trademark of Sirigen Group Ltd.</p> <p>Learn more about Brilliant Violet™.</p> <p>This product is subject to proprietary rights of Sirigen Inc. and is made and sold under license from Sirigen Inc. The purchase of this product conveys to the buyer a non-transferable right to use the purchased product for research purposes only. This product may not be resold or incorporated in any manner into another product for resale. Any use for therapeutics or diagnostics is strictly prohibited. This product is covered by U.S. Patent(s), pending patent applications and foreign equivalents.</p>
Excitation Laser	Violet Laser (405 nm)
Application Notes	Additional reported applications (for the relevant formats) include: Western blotting ¹ , immunoprecipitation ¹ , immunohistochemical staining ^{2,3} of paraffin embedded ⁷ and frozen tissue sections, ELISA ⁴ , and functional assay ⁵ . The LEAF™ or Ultra-LEAF™ purified antibody is

recommended for functional assays (contact our [custom solutions team](#)).

Application References

1. Latza U, *et al.* 1994. *Eur. J. Immunol.* 24:677. (WB, IP)
2. Durkop H, *et al.* 1995. *Brit. J. Haematol.* 91:927. (IHC)
3. Durkop H, *et al.* 1997. *Brit. J. Haematol.* 98:863. (IHC)
4. Willett B, *et al.* 2007. *J. Virol.* 81:9665. (ELISA)
5. Li M and Zhang Y. *et al.* 2005. *Cell. Mol. Immunol.* 2:467. (FA)
6. Gloviczki ML, *et al.* 2012. *Clin. J. Am. Soc. Nephrol.* 8:546. [PubMed](#)
7. Domingos PL, *et al.* 2012. *An. Bras. Dermatol.* 87:851. (IHC)

RRID

AB_2616911 (BioLegend Cat. No. 350025)

AB_2616912 (BioLegend Cat. No. 350026)

Antigen Details

Structure	50 kD, TNF receptor family
Distribution	Activated T cells
Function	Receptor for OX40 ligand, co-stimulatory signal for T cell proliferation and survival.
Ligand/Receptor	OX40L (CD252)
Cell Type	T cells, Tregs
Biology Area	Apoptosis/Tumor Suppressors/Cell Death, Cell Biology, Costimulatory Molecules, Immunology
Molecular Family	CD Molecules, Immune Checkpoint Receptors
Antigen References	<ol style="list-style-type: none">1. Smith CA, <i>et al.</i> 1994. <i>Cell.</i> 76:959.2. Chen AL, <i>et al.</i> 1999. <i>Immunity.</i> 11:689.3. Croft M. 2010. <i>Annu. Rev. Immunol.</i> 28:57.4. Ruby CE, <i>et al.</i> 2009. <i>J. Immunol.</i> 183:5079.5. Klinger M, <i>et al.</i> 2009. <i>J. Immunol.</i> 182:4581.
Gene ID	7293


Related Protocols

[Cell Surface Flow Cytometry Staining Protocol](#)


Other Formats

FITC anti-human CD134 (OX40), Purified anti-human CD134 (OX40), PE anti-human CD134 (OX40), APC anti-human CD134 (OX40), PerCP/Cyanine5.5 anti-human CD134 (OX40), PE/Cyanine7 anti-human CD134 (OX40), Brilliant Violet 421™ anti-human CD134 (OX40), Purified anti-human CD134 (OX40) (Maxpar® Ready), Alexa Fluor® 647 anti-human CD134 (OX40), PE/Dazzle™ 594 anti-human CD134 (OX40), APC/Cyanine7 anti-human CD134 (OX40), Biotin anti-human CD134 (OX40), Brilliant Violet 510™ anti-human CD134 (OX40), Brilliant Violet 605™ anti-human CD134 (OX40), APC/Fire™ 750 anti-human CD134 (OX40), Brilliant Violet 711™ anti-human CD134 (OX40), TotalSeq™-A0158 anti-human CD134 (OX40), TotalSeq™-C0158 anti-human CD134 (OX40), TotalSeq™-B0158 anti-human CD134 (OX40), TotalSeq™-D0158 anti-human CD134 (OX40), PE/Cyanine5 anti-human CD134 (OX40), PE/Fire™ 810 anti-human CD134 (OX40)

Product Data

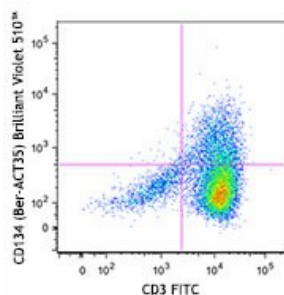

PHA-stimulated (3 days) human peripheral blood lymphocytes were stained with CD3 FITC and CD134 (OX40) (clone Ber-ACT35) Brilliant Violet 510™ (top), or mouse IgG1, κ Brilliant Violet 510™ isotype control (bottom).

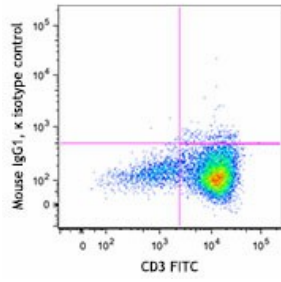

For research use only. Not for diagnostic use. Not for resale. BioLegend will not be held responsible for patent infringement or other violations that may occur with the use of our products.

*These products may be covered by one or more Limited Use Label Licenses (see the BioLegend Catalog or our website, www.biolegend.com/ordering#license). BioLegend products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products, reverse engineer functionally similar materials, or to provide a service to third parties without written approval of BioLegend. By use of these products you accept the terms and conditions of all applicable Limited Use Label Licenses. Unless otherwise indicated, these products are for research use only and are not intended for human or animal diagnostic, therapeutic or commercial use.

BioLegend Inc., 8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587